




PRODUCT ADVANTAGES

- Electromagnetic force balance sensor ensures accurate and stable results
- Multi-unit conversion: g, oz, ct, lb, GN, with power-off memory
- Adjustable zero tracking and auto-tare function
- Front leveling feet with visible level bubble
- Four-sided HD glass draft shield, 240 mm operating space
- RS232 output interface, connectable to printers and computers
- Large LED backlit display for clear and intuitive reading
- Upgradeable internal calibration

PRODUCT DETAILS

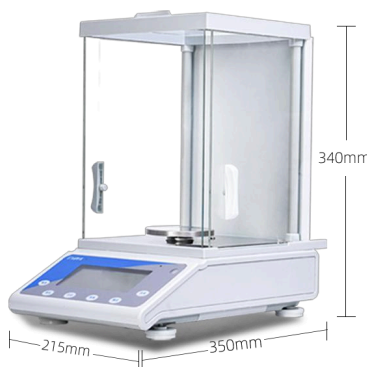

- LED DISPLAY**
 Large backlit LCD, clear reading under any lighting
- GLASS DRAFT SHIELD**
 Four-sided windshield prevents airflow interference
- LEVEL INDICATOR**
 Integrated bubble level for precise alignment
- ADJUSTABLE FEET**
 Front-accessible leveling feet for easy adjustment

PRODUCT PARAMETERS

Parameter	Standard Range	Extended Range
Weighing Range	0 - 120 g	0 - 220 g
Readability	0.1 mg	
Linearity / Repeatability	± 0.2 mg / ± 0.2 mg	
Interface / Stabilization	RS232 / ≤ 5 s	
Pan / Chamber	φ 90 mm SS / 240 mm	
Dimensions / Power	350 x 215 x 340 mm / AC 110-240V	